

www.mywam.eu



NOTE: The information in this leaflet is provided for convenience and education purposes only and therefore: a) such information does not constitute a commercial proposal (quotation) within the meaning of commercial law, b) such information cannot be the basis for taking decisions about taking up treatment, withdrawal from treatment or other medical procedures, c) the information should be interpreted by a medical professional. d) the products may be modified by the manufacturer without notice.





C € MD MEDICAL DEVICE - CLASS I

The MOUSE special needs stroller is designed for the smallest children who are unable to walk or sit in a correct position as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for children at a height of 80-110 cm and a weight of 30 kg.

Thanks to the backrest tilt and recline feature as well as adjustable torso and head pads you may easily adjust MOUSE to the required position and ensure the user with optimal protection and stability.

The shock-absorbing suspension on all four, large, gel-filled wheels together with the tilt-foot facilitate travel on uneven terrain. Modern and stylish design is an additional advantage of this stroller.

INDICATIONS

paralysis and paresis of lower limbs of various origins, paralysis and paresis of the trunk, weakening of postural muscles responsible for maintaining the correct body position, cerebral palsy, myelomening ocele, muscular dystrophies.

CONTRAINDICATIONS:

In case of any contraindications to sitting position &/or passive verticalization stroller adjustments should be consulted with physiotherapist or physician-in-charge.

colour options.

















TILT IN SPACE SEAT

simple mechanism



8°-23°

4 positions

RECLINING BACKREST

wide range of adjustments



16°-88°

4 positions

ADJUSTABLE AND ANGLE TILT FOOTREST

simple and quick adjustment



84°- 173°

6 positions

SEAT MINIMIZING KIT

reducing insert



The SEAT MINIMIZING KIT in MOUSE is a comfortable, one-peace insert that allows you to reduce both width and depth of the seat.

STABILIZATION SYSTEM

adjustable pads



The head and torso support consists of two pairs of independent pads that are height, width and angle adjustable. Its wide range of adjustments and simplicity gives you a fast and effective adaptation to the user's individual needs.

LYING POSITION

wide range of backrest recline



With a fully adjustable seat and footrest you may quickly alter the stroller into a cozy bed perfect for a short rest or a nap.

REVERSIBLE SEAT

arrange your travel easily



The design of the stroller allows you to easily rotate the seat opposite to the driving direction. Thanks to this feature you can keep an eye on your child all the time.

FOLDING LIGHTWEIGHT FRAME

ease of storage, transportation and dismantling



The construction of main elements allows you to dismantle Mouse into 3 separate units (frame, seat and wheels).

DRIVING SUPPORT

comfort for the bumpiest journeys



The stroller is equipped with a shock-absorbing suspension on all four wheels and an adjustable rear suspension in the rear wheels, swivel front castor wheels with option to lock, central brake and a tilt-foot.



accessories.wide range of equipment



parameters.

TECHNICAL DATA OF A STROLLER READY TO DRIVE



seat depth 22*/26 cm

backrest height 52 cm

footrest adjustable within 14,5 - 23 cm

- **(a)** total weight 17,2 23,2 kg
- **b** total lenght 95 124 cm
- **16** total height 122 134 cm
- total width 60 cm

dimension after folding: frame L 84 x W 58 x H 25 cm / seat L 72 x W 42 x H 39 cm

USER PARAMETERS



🐚 maximum weight 30 kg





C € MD MEDICAL DEVICE - CLASS I

The MEWA special needs stroller is designed for children who are unable to walk or sit in a correct position as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for children at a height of 90-125 cm and a weight of 40 kg.

Thanks to the backrest tilt and recline feature as well as adjustable torso and head pads you may easily adjust MEWA to the required position and ensure the user with optimal protection and stability.

The shock-absorbing suspension on all four, large, gel-filled wheels together with the tilt-foot facilitate travel on uneven terrain. Modern and stylish design is an additional advantage of this stroller.

INDICATIONS

paralysis and paresis of lower limbs of various origins, paralysis and paresis of the trunk, weakening of postural muscles responsible for maintaining the correct body position, cerebral palsy, myelomeningocele, muscular dystrophies.

CONTRAINDICATIONS:

In case of any contraindications to sitting position &/or passive verticalization stroller adjustments should be consulted with physiotherapist or physician-in-charge.

colour options.



UPHOLSTERY



FRAIVIE spray painted / glossy



UPHOLSTERY



FRAME spray painted / half matt













FRAME spray painted / glossy



UPHOLSTERY synthetic



FRAME
spray painted / gloss



UPHOLSTERYsynthetic



FRAME spray painted / glossy













TILT IN SPACE SEAT

simple mechanism



5,5°-21,5°

4 positions

RECLINING BACKREST

wide range of adjustments



8°-86°

4 positions

ADJUSTABLE AND ANGLE TILT FOOTREST

simple and quick adjustment



64°- 180°

5 positions

SEAT MINIMIZING KIT

reducing inserts



The SEAT MINIMIZING KIT in MEWA is a set of 4 comfortable, independent inserts that allows you to reduce both width and depth of the seat.

STABILIZATION SYSTEM

adjustable pads



The head and torso support consists of two pairs of independent pads that are height, width and angle adjustable. Its wide range of adjustments and simplicity gives you a fast and effective adaptation to the user's individual needs.

LYING POSITION

wide range of backrest recline



With a fully adjustable seat and footrest you may quickly alter the stroller into a cozy bed perfect for a short rest or a nap.

REVERSIBLE SEAT

arrange your travel easily



The design of the stroller allows you to easily rotate the seat opposite to the driving direction. Thanks to this feature you can keep an eye on your child all the time

FOLDING LIGHTWEIGHT FRAME

ease of storage, transportation and dismantling



The construction of main elements allows you to dismantle MEWA into 3 separate units (frame, seat and wheels).

DRIVING SUPPORT

comfort for the bumpiest journeys



The stroller is equipped with a shock-absorbing suspension on all four wheels and an adjustable rear suspension in the rear wheels, swivel front castor wheels with option to lock, central brake and a tilt-foot.



accessories.wide range of equipment



parameters.

TECHNICAL DATA OF A STROLLER READY TO DRIVE



seat depth 24* / 28 cm

backrest height 62 cm

😽 footrest adjustable within 17 - 32 cm









dimension after folding: frame L 84 x W 64 x H 26 cm / seat L 88 x W 43 x H 42 cm

USER PARAMETERS









C € MD MEDICAL DEVICE - CLASS I

The GRIZZLY special needs stroller is designed for the children who are unable to walk or sit in a correct position as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for children at a height of 85-140 cm and a weight of 45 kg.

Grizzly is equipped with the "I GROW WITH YOU" system that allows for a smooth extension of the backrest height and seat depth as your child grows. Moreover Grizzly offers you a comprehensive stabilization system consisting of adjustable pads at the level of the head, hips and chest – together with a fully adjustable footrest, tilt in space seat and reclining backrest we bring you a stroller that is tailored to fit your child's needs and will serve for years. Grizzly is crash tested and it ensures safety while transporting children in specially adapted vehicles.

INDICATIONS

paralysis and paresis of lower limbs of various origins, paralysis and paresis of the trunk, weakening of postural muscles responsible for maintaining the correct body position, cerebral palsy, myelomening ocele, muscular dystrophies.

CONTRAINDICATIONS:

In case of any contraindications to sitting position &/or passive verticalization stroller adjustments should be consulted with physiotherapist or physician-in-charge.

colour options.



UPHOLSTERY



spray painted / glossy




UPHOLSTERY synthatic

FRAME
spray painted / glossy







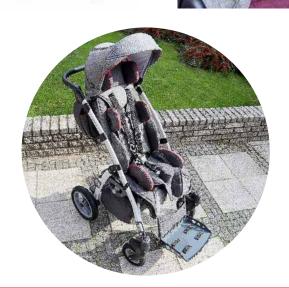

TILT IN SPACE SEAT simple mechanism



0°-30° 4 positions

I GROW WITH YOU extensive adjustment options



"I GROW WITH YOU" system allows to extend:

- the seat depth up to 37 cm
- backrest height up to 70 cm (no tools required) Designed to adapt and evolve!

RECLINING BACKREST

wide range of adjustments



91°-212° 6 positions

STABILIZATION SYSTEM

safety and stability



A complete support system consisting of a four independent pairs of pads (head, torso, hips and lateral) guarantees maximum of adjustment possibilities as well as highest comfort and safety during travel.

ADJUSTABLE AND ANGLE TILT FOOTREST

simple and quick adjustment



55°- 215° 5 positions

LYING POSITION

wide range of backrest recline



With a fully adjustable seat and footrest you may quickly alter the stroller into a cozy bed perfect for a short rest or a nap.

REVERSIBLE SEAT

arrange your travel easily



The design of the stroller allows you to easily rotate the seat opposite to the driving direction. Thanks to this feature you can keep an eye on your child all the time.

FOLDING LIGHTWEIGHT FRAME

ease of storage, transportation and dismantling



The construction of main elements allows you to dismantle GRIZZLY into 3 separate units (frame, seat and wheels).

DRIVING SUPPORT

comfort for the bumpiest journeys



The stroller is equipped with a shock-absorbing suspension on all four wheels and an adjustable rear suspension in the rear wheels, swivel front castor wheels with option to lock, central brake and a tilt-foot.



accessories. wide range of equipment



parameters.

TECHNICAL DATA OF A STROLLER READY TO DRIVE



seat depth 25 - 37 cm

backrest height 57 - 70 cm

🈽 footrest adjustable within 23 - 40,5 cm



total lenght 82 - 135 cm

🔈 total height 108 - 125 cm


total width 65 cm

dimension after folding: frame L96 x W 65 x H 51 cm / seat L72 x W 43 x H 47 cm

USER PARAMETERS



mayimum waight 15

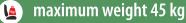





C € MD MEDICAL DEVICE - CLASS I

The Yeti special needs stroller is designed for larger/older children who are unable to walk or sit in a correct position as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for users at a height of 130-175 cm and a weight of 80 kg.

Its heavy-duty construction is equipped with the "I GROW WITH YOU" system (smooth extension of the backrest height and the seat depth as your child grows) and adjustable tension system which allows to adjust the proper upholstery tension both on the seat and the backrest guarantees you a support for many years.

Moreover, Yeti's solid suspension along with a fast folding mechanism, tilted seat, reclining backrest and a wide range of accessories significantly simplify the everyday use irrespective of the weather conditions and time of the year. Modern and practical design is his additional asset.

INDICATIONS

paralysis and paresis of lower limbs of various origins, paralysis and paresis of the trunk, weakening of postural muscles responsible for maintaining the correct body position, cerebral palsy, myelomening ocele, muscular dystrophies.

CONTRAINDICATIONS:

In case of any contraindications to sitting position &/or passive verticalization stroller adjustments should be consulted with physiotherapist or physician-in-charge.

colour options.





UPHOLSTERY







UPHOLSTERY



spray painted / matt















TILT IN SPACE SEAT

simple mechanism

RECLINING BACKREST

wide range of adjustments

ADJUSTABLE AND ANGLE TILT FOOTREST

simple and quick adjustment







9,5°- 9°

5 positions

3°- 99°

8 positions

80°-180°

6 positions

SEAT ADJUSTMENTS

reducer cushions



Reducer cushions set consisting of 4 independent velcro pads offers an additional adaptation of the seat interior through simple attaching of the pads at desired position. Easy adjustment of seat width and depth at your fingertips.

I GROW WITH YOU

extensive adjustment options



- "I GROW WITH YOU" system allows to extend:
- the seat depth up to 45 cm
- backrest height up to 85 cm (tools required) Designed to adapt and evolve!

LYING POSITION

wide range of backrest recline



With a fully adjustable seat and footrest you may quickly alter the stroller into a cozy bed perfect for a short rest or a nap.

WINDSTOPPER

multifunctional upholstery



YETI is equipped with a multifunction windstopper that not only comes in handy during windy conditions but also protects against unnecessary arm movement and limits visual stimulation.

FOLDING LIGHTWEIGHT FRAME

ease of storage and transportation



The one-piece construction of YETI stroller allows for a quick and easy folding which facilitates daily use, storage and transport.

DRIVING SUPPORT

comfort for the bumpiest journeys



The stroller is equipped with adjustable rear suspension in the rear wheels, swivel front castor wheels with option to lock, central brake and a tilt-foot.



accessories.wide range of equipment



parameters.

TECHNICAL DATA OF A STROLLER READY TO DRIVE



seat depth 31 - 45 cm

backrest height 65 - 85 cm

🈽 footrest adjustable within 28,5 - 43,5 cm 👔











USER PARAMETERS









C € MD MEDICAL DEVICE - CLASS I

The PEGAZ special needs stroller is designed for the children who are unable to walk or sit in a correct position as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for users at a height of 80-140 cm and a weight of 40 kg.

PEGAZ is an ultralight umbrella type stroller made of aluminium and steel alloys equipped with the "I GROW WITH YOU" system, narrowing inserts, reclining backrest and a wide range of accessories. Simple to fold and unfold, conveniently compact (staring from 10 kg) and easy to manoeuvre – Pegaz will be your day-to-day support for many years.

paralysis and paresis of lower limbs of various origins paralysis and paresis of the trunk weakening of postural muscles responsible for maintaining the correct body position cerebral palsy myelomeningocele muscular dystrophies

Contraindications:

contraindications to sitting position contraindications to passive verticalization

colour options.

























SEAT POSITION

angle parameter



25°

1 position

RECLINING BACKREST

angles parameters



22°, 34°

2 positions

ADJUSTABLE FOOTREST

hight regulation



20, 30, 40 cm

3 positions

I GROW WITH YOU

extensive adjustment options



- "I GROW WITH YOU" system allows to:
- extend the backrest height by 10 cm, giving you maximum height of 60 cm $\,$
- extend the seat depth by 8 cm, giving you maximum depth of 33 cm.

 Quick and simple adjustment designed to adapt and evolve!

SEAT ADJUSTMENTS

narrowing inserts



PEGAZ is equipped with 2 narrowing inserts that allow to reduce the seat width and provide additional safety and stability during travel.

FOLDING LIGHTWEIGHT FRAME

ease of storage and transportation

DRIVING SUPPORT

comfort and maneuverability









The one-piece, umbrella construction of PEGAZ stroller allows for a quick and easy folding which due to its compact size and low weight facilitates daily use, storage and transport.

The stroller is equipped with a shock-absorbing suspension on all four wheels, swivel front castor wheels with option to lock, two independent rear wheel brakes and a tilt-foot.



accessories.wide range of equipment



parameters.


TECHNICAL DATA OF A STROLLER READY TO DRIVE



seat depth 25/33 cm

backrest height 50/60 cm

footrest adjustable within 20/30/40 cm

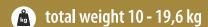








USER PARAMETERS



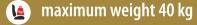









C € MD MEDICAL DEVICE - CLASS I

The ZEFIR special needs stroller is designed for the children who are unable to walk or sit in a correct position as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for users at a height of 80-140 cm and a weight of 40 kg.

ZEFIR is a lightweight stroller made of aluminium and steel alloys equipped with the "I GROW WITH YOU" system, narrowing inserts, reclining backrest and a wide range of additional accessories. Simple to fold and unfold, conveniently compact and easy to manoeuvre – ZEFIR will be your dayto-day support for many years.

INDICATIONS:

paralysis and paresis of lower limbs of various origins paralysis and paresis of the trunk weakening of postural muscles responsible for maintaining the correct body position cerebral palsy myelomeningocele muscular dystrophies

CONTRAINDICATIONS:

contraindications to sitting position contraindications to passive verticalization

colour options.



UPHOLSTERY



FRAME





UPHOLSTERY





FRAME

spray painted / matt

















SEAT POSITION

angle parameter



25°

1 position

RECLINING BACKREST

angles parameters



25°, 57°

2 positions

ADJUSTABLE FOOTREST

hight regulation



20, 30, 40 cm

3 positions

I GROW WITH YOU

extensive adjustment options



- "I GROW WITH YOU" system allows to:
- extend the backrest height by 10 cm, giving you maximum height of 60 cm $\,$
- extend the seat depth by 6 cm, giving you maximum depth of 30 cm.

 Quick and simple adjustment designed to adapt and evolve!

SEAT ADJUSTMENTS

narrowing inserts



ZEFIR is equipped with 2 narrowing inserts that allow to reduce the seat width and provide additional safety and stability during travel.

FOLDING LIGHTWEIGHT FRAME

ease of storage and transportation

DRIVING SUPPORT

comfort and maneuverability



Light and compact construction of the Zefir special needs stroller allows you to quickly pack it and take it anywhere you want! Enjoy easy manoeuvre across both the city and scenic all-terrain routes.





The special needs stroller Zefir is equipped with front wheels swivel lock, shock absorption system on the rear wheels and central brake mounted on the rear beam. Next to the brake we have mounted a tilt-foot which helps to lift front of the stroller while driving through curbs and other small obstacles. One of Zefir's great advantage is the possibility to mount the all-terrain wheels that allow you to explore new places and enjoy outdoor activities like never before.



akcesoria. bogate możliwości wyposażenia



parameters.

TECHNICAL DATA OF A STROLLER READY TO DRIVE

- seat width 29*/38 cm
- seat depth 24/30 cm
- backrest height 50/60 cm
- b⇒ footrest adjustable 20/30/40 cm
- kg
 - total weight 13 16 / ** 13 18 kg
 - 📐 total lenght 96 111 / ** 98 -113 cm
- 🔈 total height 100 136 / ** 102 138 cm
- total width 60 / ** 62 cm

USER PARAMETERS

- - growth 80 140 cm
- maximum weight 40 kg
 - * * with reducer cushions
 - ** with all-terrain wheels



dimension after folding: L 93 x W 59 x H 26 cm / **L 85 x W 60 x H 31 cm



www.mywam.eu

contact:

EN mob. +48 530 304 144 partner@mywam.pl

ES / EN mob. + 48 530 022 520 export@mywam.pl RUS / EN mob. + 48 575 997 180 yuliia@mywam.pl